

DRUG-FREE **BACK PAIN RELIEF**



WiTouch Pro
BLUETOOTH® TENS THERAPY
KEEP MOVING™



 **Bluetooth®**

 Download on the
App Store

 GET IT ON
Google Play

- 
- Deep penetrating pain relief right where you need it
 - Doctor office quality
 - Blocks the body's pain signal and releases endorphins
 - Alternative to pharmaceutical opioids



Tens Therapy



Made in the US




FDA approved

Wi Touch Pro uses Transcutaneous Electrical Nerve Stimulation (TENS) to diminish pain and replace it with a tingling or massage-like sensation. It's a safe, simple and effective therapy that can be applied to the upper, middle or lower back using gel pad adhesives.

■ Specifications:

- Waveform
- Asymmetrical Biphasic
- Output: 55V Max
- 120Hz Max Frequency
- 250Qs Max pulse Duration

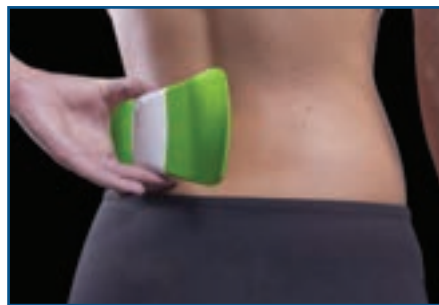

Program	Description	Frequency	Pulse Duration
WiTouch Exclusive	High/Low frequency combination for benefits of both opiate Release and Gate Control theories for TENS pain relief	80-120Hz 5-10Hz	120-μsec 240-μsec
High Frequency	High frequency only to provide benefit of Gate control, providing a localized, soft flowing stimulation sensation	80-120Hz	120-μsec
Low Frequency	Low frequency only to provide benefit of Opiate Release, providing a pulsing stimulation sensation	50-10Hz	240-μsec
High-Low Combo	Low frequency alternating with high frequency to provide benefits of both Opiate Release and Gate Control	2Hz-70Hz	250-μsec 150-μsec